




CLEANTYPE® PRIME SERIES

TOUCH+ (WITH BACKLIGHT)
TOUCH

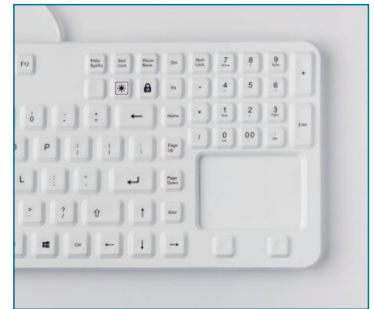

SPECIAL FEATURES

- IP68 rated medical keyboard
- Adjustable brightness of backlight (TOUCH+)
- Fully washable and disinfectable
- 109 keys incl. num pad and track pad
- Waterproof silicone USB protection cap
- Excellent tactile feeling
- Integrated Clean Mode for cleaning
- Magnetic mounting
- EN60601-1-2 Ed 4 certified

SPECIFICATION

Number of keys:	109 (US-Layout)	Interfaces:	USB 2.0
Keyswitch technology:	carbon contact	Cable length:	15cm
Switching cycles:	approx. 2 mill. (per key)	Operating temperature:	0°C ~ +55°C
Protection level:	IP68	Storage temperature:	-30°C ~ +70°C
Housing material:	silicone	Dimensions:	365(L) x 134(W) x 14.2(T) mm
Housing color:	white & black	Weight:	approx. 615g
Backlight:	blue color, 5 level adjustable brightness	Certification:	EN60601-1-2 Ed 4, CE, FCC, RoHS, REACH
		Warranty:	2 years

ORDER INFORMATION

Cat. No.	Product description	Colour	Layout
KSI-U10210	CLEANTYPE® PRIME TOUCH+	White	US/DE/JP/KR/THAI
KSI-U10211	CLEANTYPE® PRIME TOUCH	White	US/DE/JP/KR/THAI
KSI-U20210	CLEANTYPE® PRIME TOUCH+	Black	US/DE/JP/KR/THAI
KSI-U20211	CLEANTYPE® PRIME TOUCH	Black	US/DE/JP/KR/THAI

Other layouts, configurations and interfaces on request.
Customized key function on request (Firmware update via USB)



EN60601-1-2 Ed 4 certified

EN 60601-1-2 Ed 4 applies to the basic safety and essential performance of Medical Equipment (ME) equipment and ME systems in the presence of electromagnetic disturbances and to electromagnetic disturbances emitted by ME equipment and ME systems.

SHORT OPERATING INSTRUCTIONS



Press and hold lock key for 5s to enable Clean Mode. Clean Mode is indicated by breathing backlight.



Adjust backlight level with this button. Five level can be adjusted including off.

CLEANING INSTRUCTION

For further details please refer to below link:
<http://www.gett-group.com/content/service-cleaning-instructions>